

# The Profiler

Now on the **FONIX® FP40**


The Profiler, a comprehensive, automated test sequence, is available on the FONIX FP40. This test sequence is a **powerful diagnostic tool**, combining ten tests chosen for their usefulness in diagnosing hearing aid problems. It is also **fast**—start-to-finish in less than 45 seconds.

## The Profiler includes the following measurements...

- OSPL90
- Maximum output
- Total harmonic distortion
- Equivalent input noise (EIN)
- Battery current drain
- 3 use gain curves at soft, medium, and loud speech levels
- Curve 2 overall gain
- Reserve gain

## The Profiler is ideal for...

- Diagnosing aids that need to be repaired or reprogrammed.
- Identifying mystery aids.
- Providing third party payor documentation.
- Identifying aids that need to be replaced.
- Establishing baseline performance for digital and programmable aids.

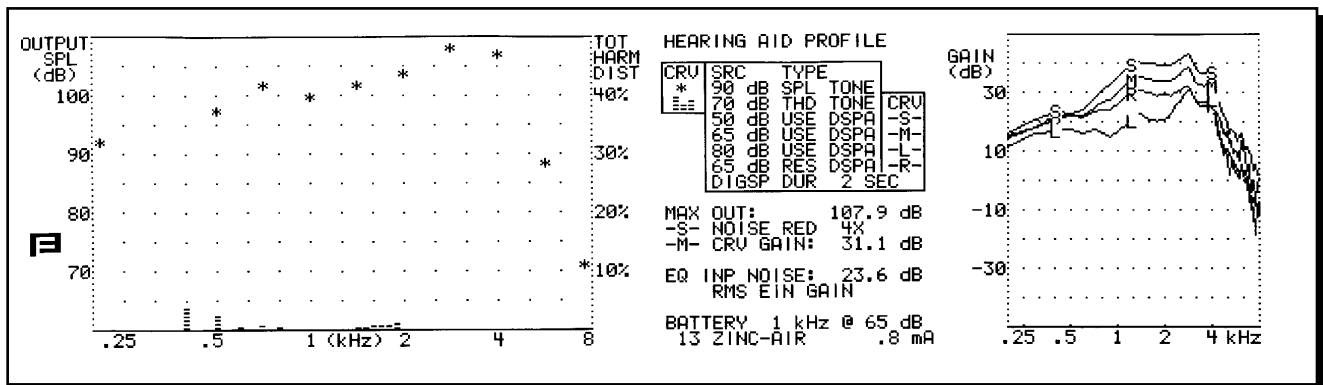

Profile of a properly functioning digital aid



**FRYE ELECTRONICS, INC.**  
P.O. Box 23391 • Tigard, OR 97281-3391  
(503) 620-2722 • (800) 547-8209  
Fax: (503) 639-0128  
www.frye.com • e-mail: sales@frye.com