


The Dual Calibration Option

for the FONIX® FA-10 and FA-12

- With the Dual Calibration Option, the FA-10/FA-12 will support two calibrations: one for TDH39 headphones and one for 3A insert earphones.

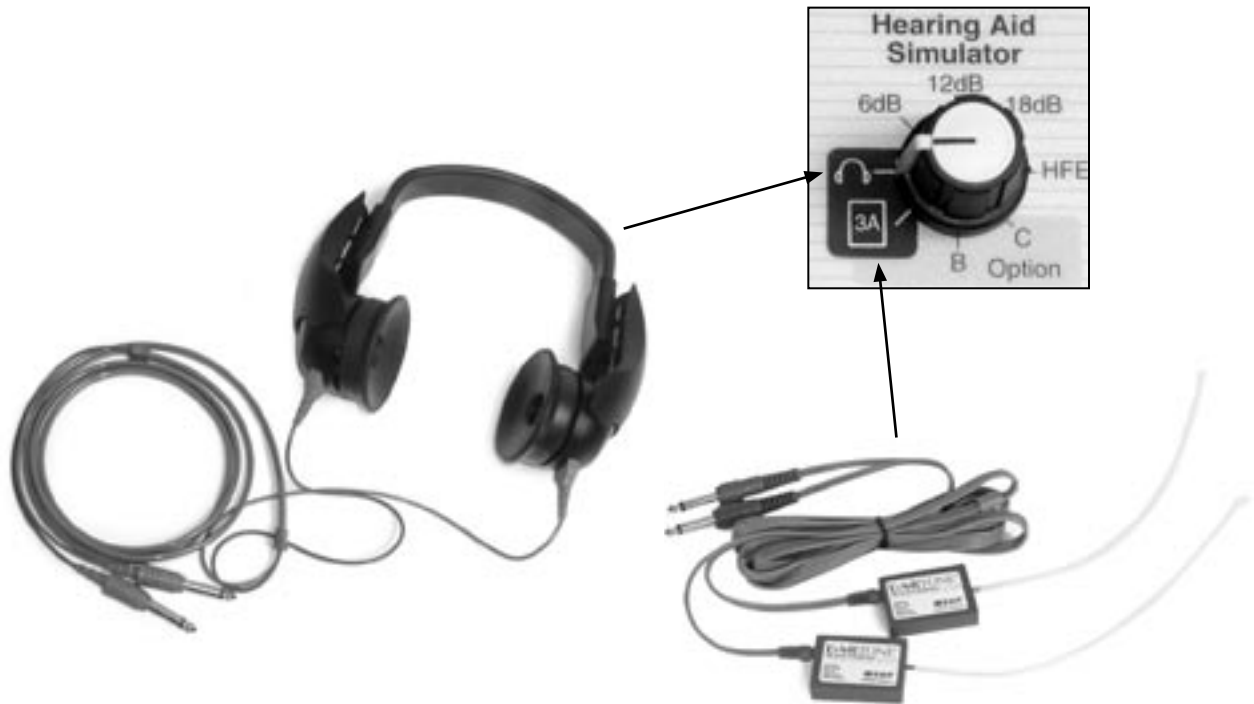

- Upon purchasing the audiometer, the customer selects a default calibration. This can be either the TDH-39 or the insert earphones. The other calibration is considered the alternate calibration.
- When the alternate calibration is selected, the right hearing aid simulator knob controls both audiometer channels. Otherwise, the right and left hearing aid simulator knobs can be controlled independently.
- The Option A, B, and C tests (ABLB, MLB, SISI) can be run using either calibration.



FRYE ELECTRONICS, INC.
P.O. Box 23391 • Tigard, OR 97281-3391 • USA
(503) 620-2722 • (800) 547-8209
Fax: (503) 639-0128
www.frye.com • e-mail: sales@frye.com

REF: DOCUMENTS/ BROCHURES/ FA-10/ DUALCAL
REV: 8/29/06