


FONIX 8000 VISIBLE SPEECH

Now with Built-in Test Signals!

We have recently updated the FONIX 8000 Visible Speech signal to include built-in real speech signals! These signals are scientifically calibrated during the leveling process in order to create accurate, repeatable measurements.

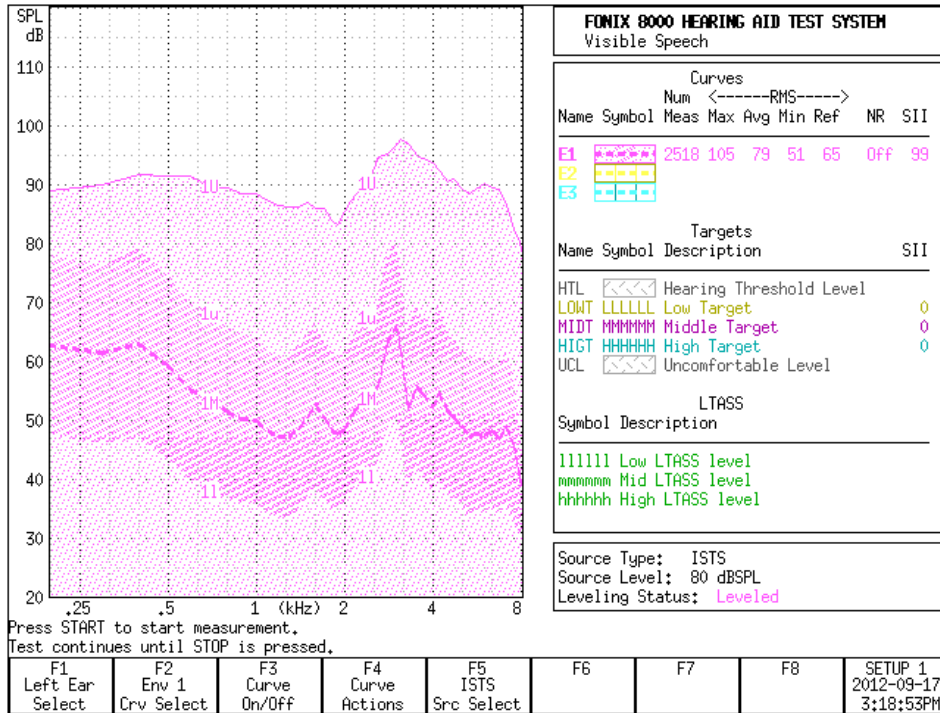

Several different signals are available including Male speaker, Female speaker, ISTS (International Speech Test Signal), White Noise, and Pink Noise. Live Voice is also available by using the External signal selection.

The FONIX 8000 Visible Speech test screen gives the clinician a complete picture of the real-ear response of the hearing aid to a built-in speech signal or to live voice. With Visible Speech, you can see:

- The patient's audiometric information (HTLs and UCLs)
- Real-ear targets at soft, medium, and loud levels (NAL-NL1, MOD-NAL)
- Real-time response to the speech signal
- Average response to the speech signal
- Speech envelope around the average response
- Minimum and Maximum response (optional)
- Reference microphone measurement (optional)
- Speech Intelligibility Index (SII)

This is a free update for all FONIX 8000 Hearing Aid Test Systems with the Real-ear Option. Just email us with your serial number and we will send you an update.



FRYE ELECTRONICS, INC.
 P.O. Box 23391 • Tigard, OR 97281-3391 • USA
 (503) 620-2722 • (800) 547-8209
 Fax: (503) 639-0128
 www.frye.com • email: sales@frye.com