


Exterior Terminology and Graphics HEARLab System

Release Date Document Change Log


September 15, 2009 Initial Release

Controlled Document - CONFIDENTIAL Effective July 28, 2009
\\sauron\engineering\design development\HEARLab\Technical File\Design Outputs\HEARLab Exterior
Terminology and Graphics.pdf page 1 of 1

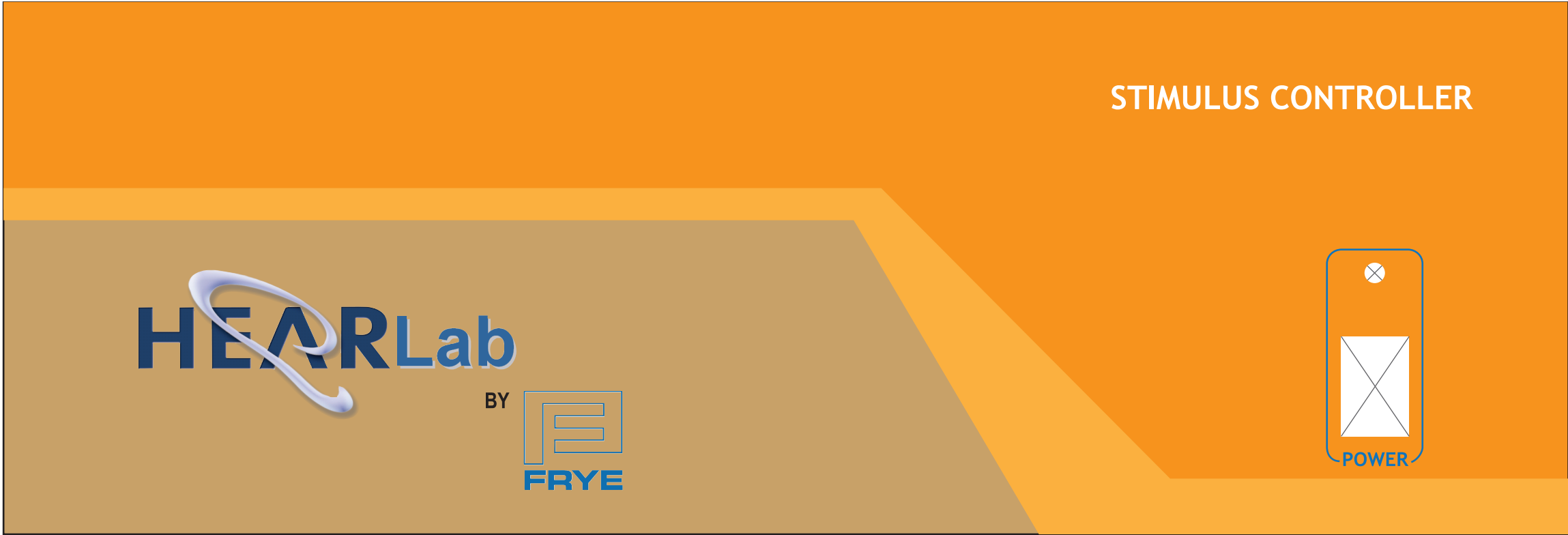

Color: CMYK
Material: GE Lexan 8B35
(Velvet Texture)
3M 9460 VHB

FRYE ELECTRONICS, INC. 9826 S.W. TIGARD ST. TIGARD, OREGON 97223	
TITLE:	Panel, Front, HEARLab
DRAWING NO:	999-9567-00
PART NO.	
DATE:	6/20/08

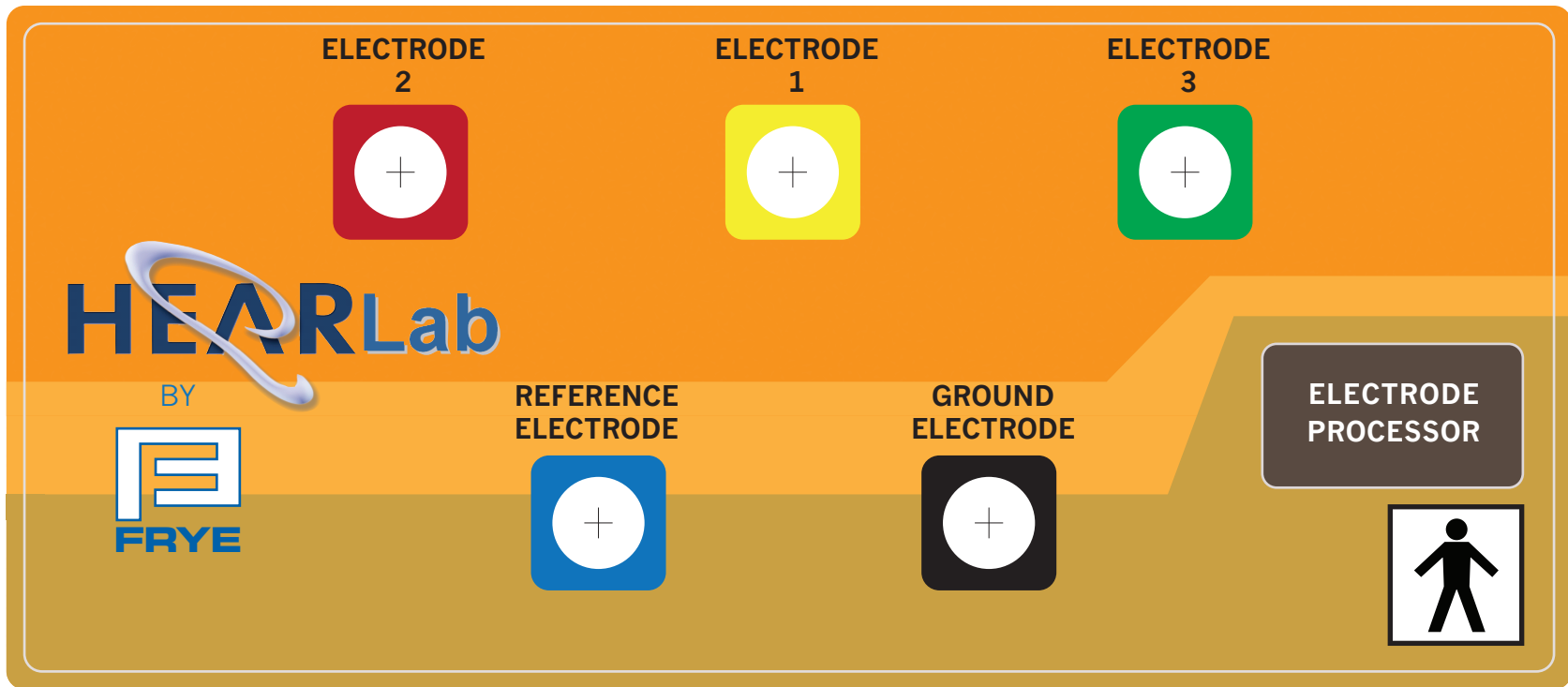

Size: 8.7 in. x 3.8 in.

Color: CMYK

Material: KrockLock+Lexan

Die cut: 5 - ø13mm

FRYE ELECTRONICS, INC.
9826 S.W. TIGARD ST. TIGARD, OREGON 97223

TITLE: Panel, Front, EP, HEARLab

DRAWING NO: 999-9568-01

PART NO. 093-0386-01

DATE: 10/22/08

Distractor Headphones

Line out

Insert Right

Headphone Right

Speaker Right

To Amplifier

From Amplifier

Speaker Left

Headphone Left

Bone

Insert Left


Monitor Speaker

+

ATTENTION: Ground reliability can be achieved ONLY when equipment is connected to “HOSPITAL GRADE” receptacles. For service refer to qualified personnel only. Consult accompanying documents.



+

 Made in Tigard, Oregon, USA
By Frye Electronics, Inc.
HEARLab is a registered trademark of HearWorks Pty Ltd
Patents Pending.

+

Caution
Amplifier outputs elevated 6VDC above chassis ground
see instruction manual.

ELECTRICAL RATING

100-240 V ~

50-60 Hz .5A

To Electrode Processor

USB from computer

External Control

Auxiliary in

Response in

Control Microphone

Line in

Speech Microphone



+

STIMULUS CONTROLLER

No user-serviceable parts inside.
Do not open enclosure.

Material: Silkscreen on metal


FRYE ELECTRONICS, INC. 9826 S.W. TIGARD ST. TIGARD, OREGON 97223	
TITLE:	Panel, Rear, HEARLab
DRAWING NO:	999-9566-03
PART NO.	034-0952-03
DATE:	10/22/09

Manufactured in 2009

FRYE ELECTRONICS, INC. 9826 S.W. TIGARD ST. TIGARD, OREGON 97223	
LABEL TITLE: Manufactured in 2009	
DRAWING NO: 999-9052-0X	
PART NO.	093-0150-11
DATE:	09/02/09




FRYE ELECTRONICS, INC.	
9826 S.W. TIGARD ST. TIGARD, OREGON 97223	
LABEL TITLE:	SAFETY GND
DRAWING NO:	999-9021-00
PART NO.	093-0222-00

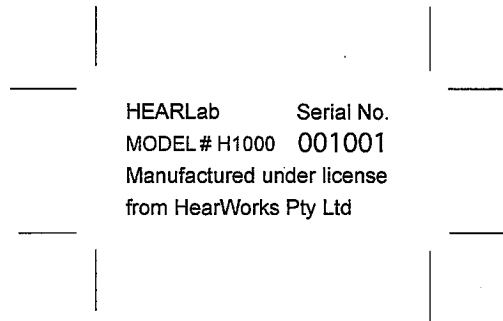


Size: 1.75 in. x 0.75 in. Round Cornered Rectangle
The incrementing serial number. Start with 001001
Material: white vinyl w/matte laminate

FRYE ELECTRONICS, INC. 9826 S.W. TIGARD ST. TIGARD, OREGON 97223	
TITLE:	Serial No., Panel, Rear, HEARLab
DRAWING NO:	999-9123-00
PART NO.	093-0400-00
DATE:	3/26/09

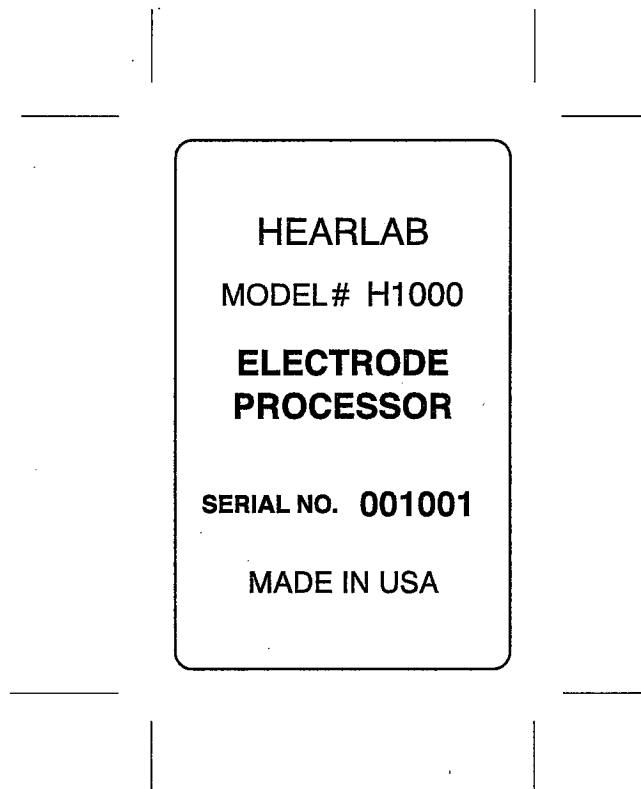

Size: 2 in. x 3 in. Round Cornered Rectangle
The incrementing serial number. Start with 001001
Material: KrockLock+Lexan

FRYE ELECTRONICS, INC. 9826 S.W. TIGARD ST. TIGARD, OREGON 97223	
TITLE:	Label, Serial No., EP, HEARLab
DRAWING NO:	999-9121-00
PART NO.	093-0397-00
DATE:	12/1/08

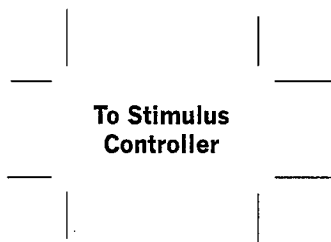

Size: 1 in. x 0.5 in. Round Cornered Rectangle
Material: KrockLock+Lexan

FRYE ELECTRONICS, INC. 9826 S.W. TIGARD ST. TIGARD, OREGON 97223	
TITLE:	Label, Side., EP, HEARLab
DRAWING NO:	999-9122-00
PART NO.	093-0399-00
DATE:	6/17/08



Kimdura, white Perm Adhesive
Size: ø 0.45 in. Circle
Color: Rich Black only
Lam: Clear laminate

FRYE ELECTRONICS, INC.
9826 S.W. TIGARD ST. TIGARD, OREGON 97223

LABEL TITLE: Electrode Earphone Frye logo

DRAWING NO: 999-9116-00

PART NO: 093-0394-00

DATE: 07/30/08



FRYE ELECTRONICS, INC. 9826 S.W. TIGARD ST. TIGARD, OREGON 97223	
LABEL TITLE:	CE
DRAWING NO:	999-9039-01
PART NO.	093-0104-01
DATE:	6/24/98

Calibrated by _____ Date Cal _____ Date Due _____

Default Tranducers

Alternate Tranducers

Right (Red) # _____ Right (Red) # _____

Left (Blue) # _____ Left (Blue) # _____

B 71 Bone Vibrator # _____

IMPORTANT

This instrument is calibrated for use with the above S/N referenced accessories.

Frye Electronics, Inc. • www.frye.com • 800-547-8209 • 503-620-2722


FRYE ELECTRONICS, INC.
9826 S.W. TIGARD ST. TIGARD, OREGON 97223

LABEL TITLE: Audiom. Calibration label

DRAWING NO: 999-9098-01


PART NO. 093-0363-01

DATE: Aug. 24, 2006

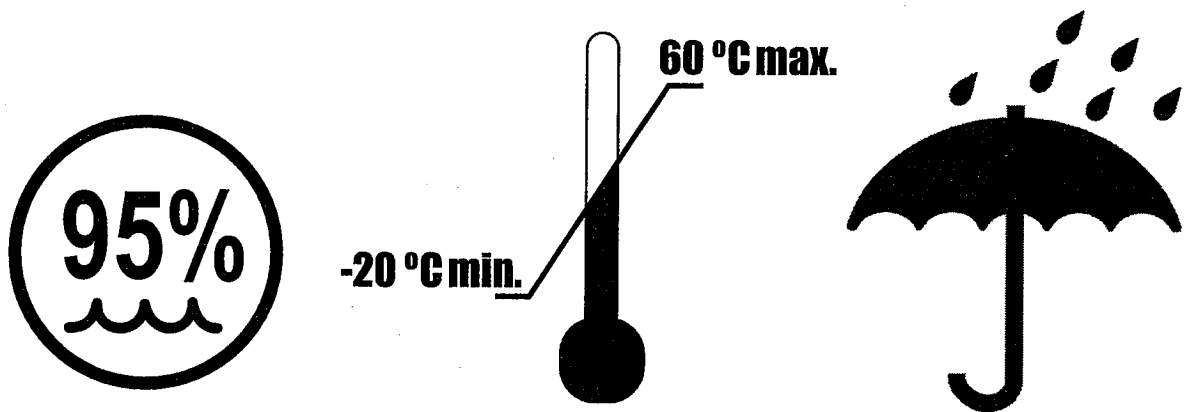

Color: Rich Black only
print on the box

FRYE ELECTRONICS, INC. 9826 S.W. TIGARD ST. TIGARD, OREGON 97223	
LABEL TITLE: Shipping Box Labeling	
DRAWING NO: 999-9117-00	
PART NO:	
DATE:	08/22/08

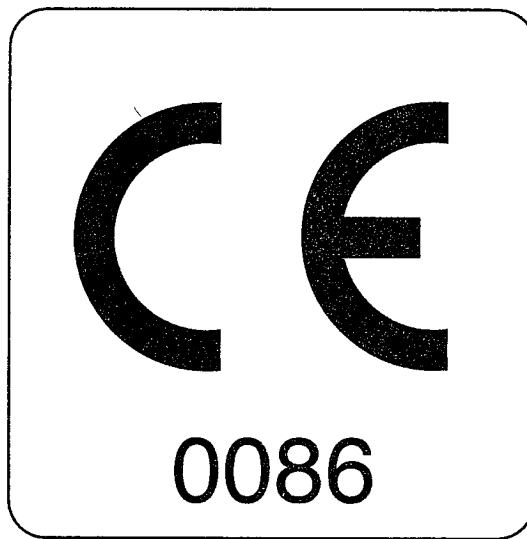

FRYE ELECTRONICS, INC.
9826 S.W. TIGARD ST. TIGARD, OREGON 97223

LABEL TITLE: CE box label

DRAWING NO: 999-9043-00

PART NO. 093-0294-00

DATE: 8/12/94